


MICROTM


INVASIVE TECHNOLOGIES
by QualiMed[®]

GWQ 18 0.018" Guidewire

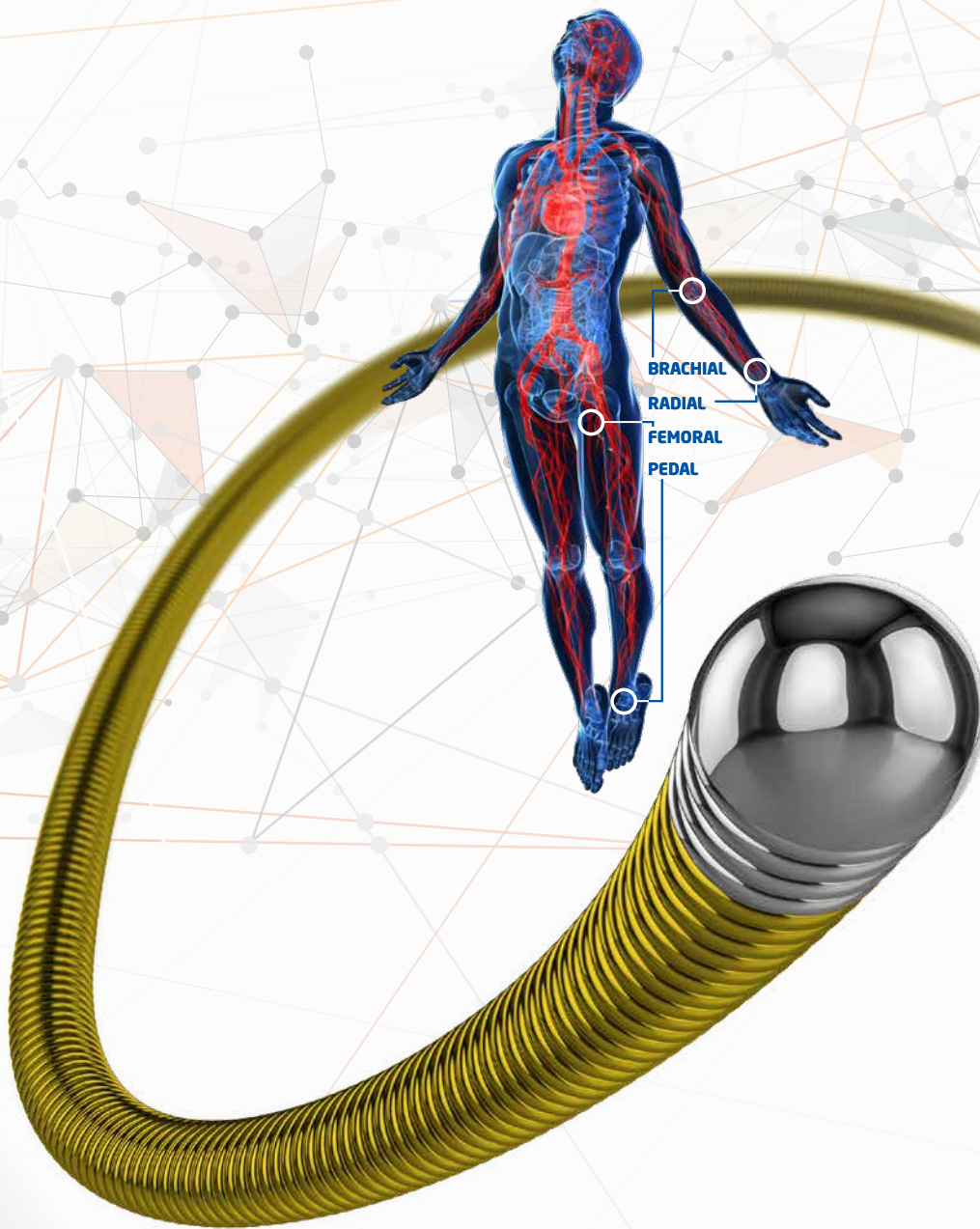

GWQ 18

0.018" Guidewire

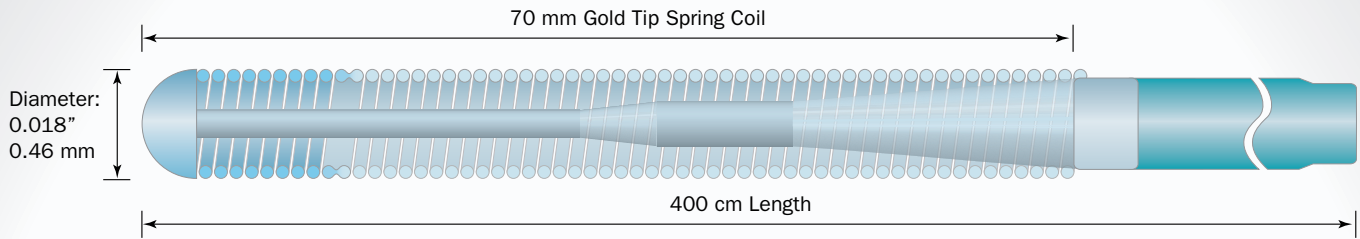

KEY FEATURES & BENEFITS

- > Mandrel guidewire construction uses a **tapered core wire** which is soldered to a coil at the distal end.
- > **Transitionless guidewire.**
- > **Gold** radiopaque tip.
- > **1:1 Torque and steerability Control** to provide precision with every turn.
- > The 0.018" guidewire line features solid stainless steel cores which provides **flexibility and kink resistance**, coiled tips for enhanced visualisation and a coating designed for maximum control and durability throughout a procedure.

TECHNICAL SPECIFICATIONS

Product name	Guidewire (GWQ)
Usable length	400 cm
Diameter	0.018" (0.46 mm)
Length gold tip distal spring coil	70 mm
Diameter distal coil	0.0175"
Core material	Stainless steel
Distal coil	Single taper stainless steel Mandrel wires
Proximal coil material	Stainless steel
Sterilization	Sterilized with ethylene oxide gas

ORDER INFORMATION

Product code	Length (cm)	Overall Length 400 cm	Diameter
GWQ018400	400		0.018"

Packaging: Dispenser and pouches; 5 pieces each per carton, sterile.